

OC-Auto Sampling Bottle 4



- For Quantitative Fecal Immunochemical Test (FIT)
- For Quantitative Fecal Calprotectin Test (FCa)

Ease, Accuracy and Hygiene

all in one compact body

Flat Surface

Allows for easy writing of patient information

Sample Buffer

Stabilizes the sample under specified conditions

Dual Aluminum Seals

Prevents accidental leakage to ensure improved hygiene

Integrated Filter

Removes fecal matter from sample solution through the squeeze filtration process

Scraper

Helps ensure accurate quantitative results by removing excess feces in the chamber

Probe Groove

Provides smooth and precise sampling of approximately 10 mg feces

Hollow cap

Easy to open and close, preventing sample leakage and contamination hazards



Original squeeze filtration process

Only filtered sample solution is dispensed for analysis.

This mechanism minimizes carry-over and analyzer maintenance, as well as helps achieve high reproducibility.

OC-Auto Sampling Bottle 4

For Fecal Collection and Sample Return



Collection Paper

- Printed sampling instruction
- Flushable paper
- Available in different designs




Special Envelope


- Integrated absorbent pad
- Available for both 1 and 2 day sampling method

Stability in the Buffer




Hemoglobin

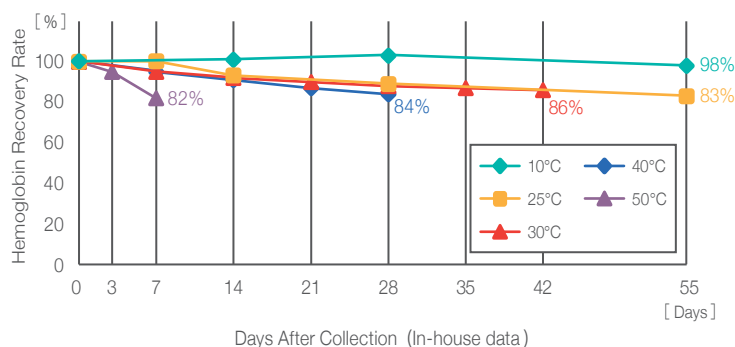

Calprotectin

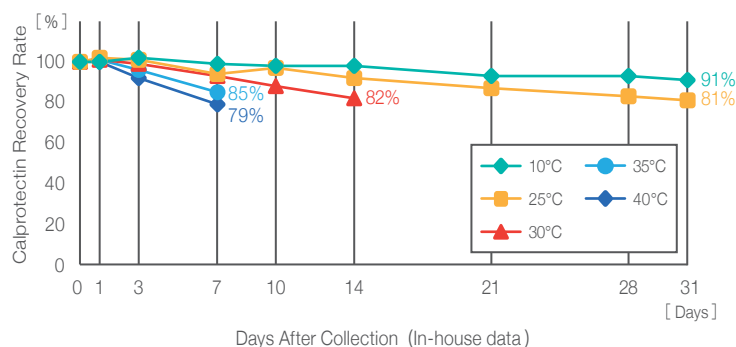

How to Collect Sample

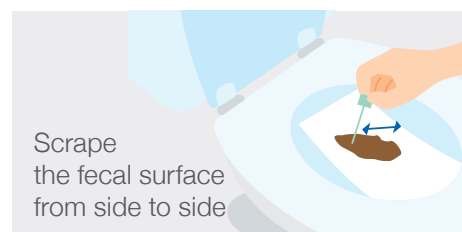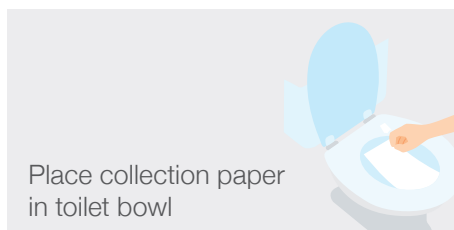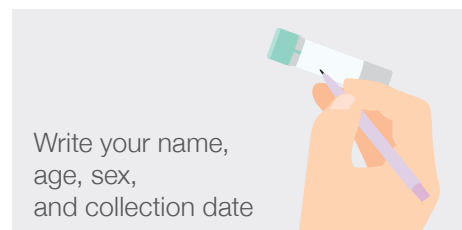

Ref No.	Product name	Contents	Storage
VCPZ45	OC-Auto Sampling Bottle 4	100 bottles	1-30°C
VCPZ46	OC-Auto Sampling Bottle 4 without barcode	100 bottles	1-30°C

Also available in other packaging options. Please contact our authorized distributor for details.

For information on use warnings and precautions, please refer to the Instructions for Use.
For Healthcare Professionals Only. Not for Public Distribution.